


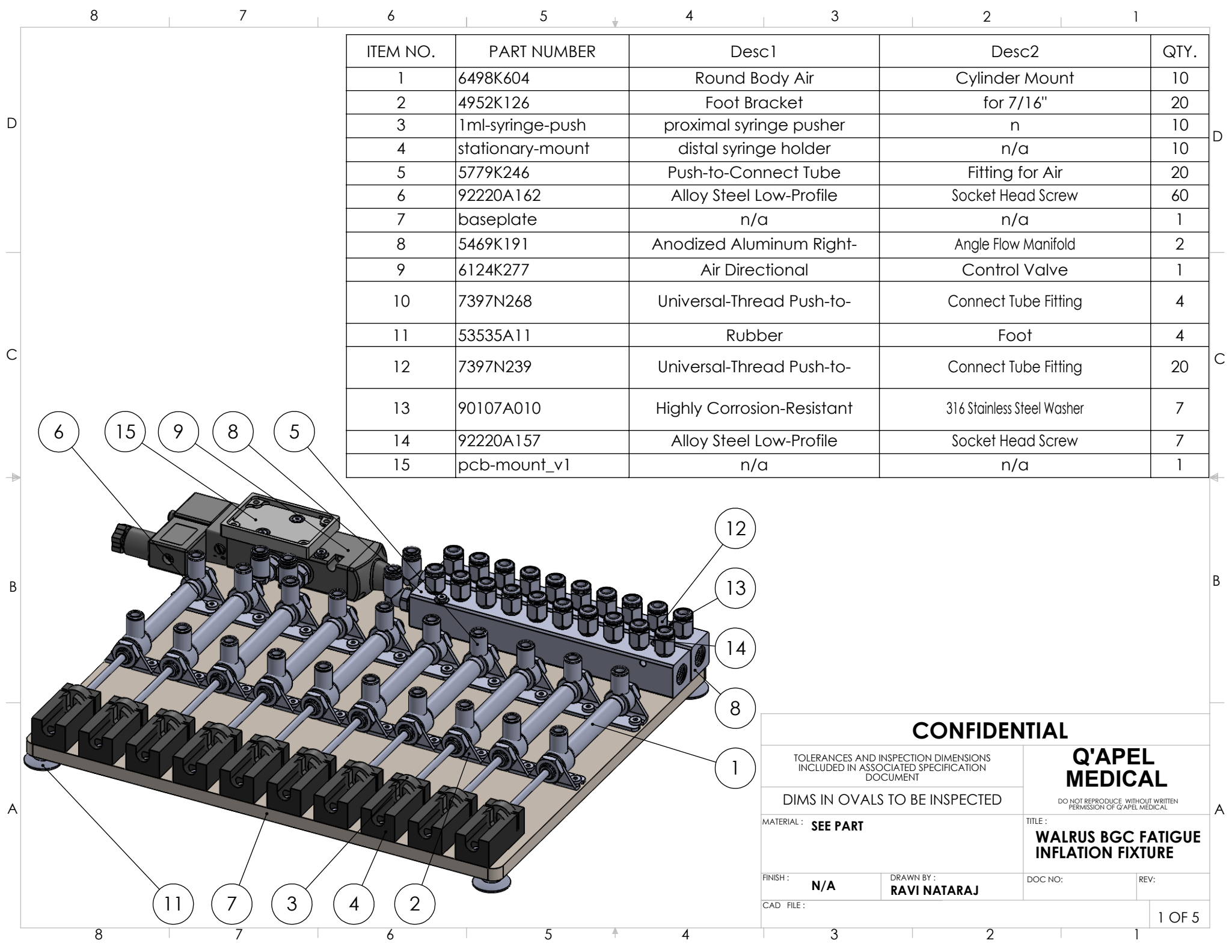

ITEM NO.	PART NUMBER	Desc1	Desc2	QTY.
1	6498K604	Round Body Air	Cylinder Mount	10
2	4952K126	Foot Bracket	for 7/16"	20
3	1ml-syringe-push	proximal syringe pusher	n	10
4	stationary-mount	distal syringe holder	n/a	10
5	5779K246	Push-to-Connect Tube	Fitting for Air	20
6	92220A162	Alloy Steel Low-Profile	Socket Head Screw	60
7	baseplate	n/a	n/a	1
8	5469K191	Anodized Aluminum Right-	Angle Flow Manifold	2
9	6124K277	Air Directional	Control Valve	1
10	7397N268	Universal-Thread Push-to-	Connect Tube Fitting	4
11	53535A11	Rubber	Foot	4
12	7397N239	Universal-Thread Push-to-	Connect Tube Fitting	20
13	90107A010	Highly Corrosion-Resistant	316 Stainless Steel Washer	7
14	92220A157	Alloy Steel Low-Profile	Socket Head Screw	7
15	pcb-mount_v1	n/a	n/a	1

CONFIDENTIAL

TOLERANCES AND INSPECTION DIMENSIONS INCLUDED IN ASSOCIATED SPECIFICATION DOCUMENT

DIMS IN OVALS TO BE INSPECTED

MATERIAL : **SEE PART**

FINISH : **N/A**

DRAWN BY : **RAVI NATARAJ**

CAD FILE :


Q'APEL MEDICAL

DO NOT REPRODUCE WITHOUT WRITTEN PERMISSION OF Q'APEL MEDICAL

TITLE : **WALRUS BGC FATIGUE INFLATION FIXTURE**

DOC NO: REV:



1 OF 5

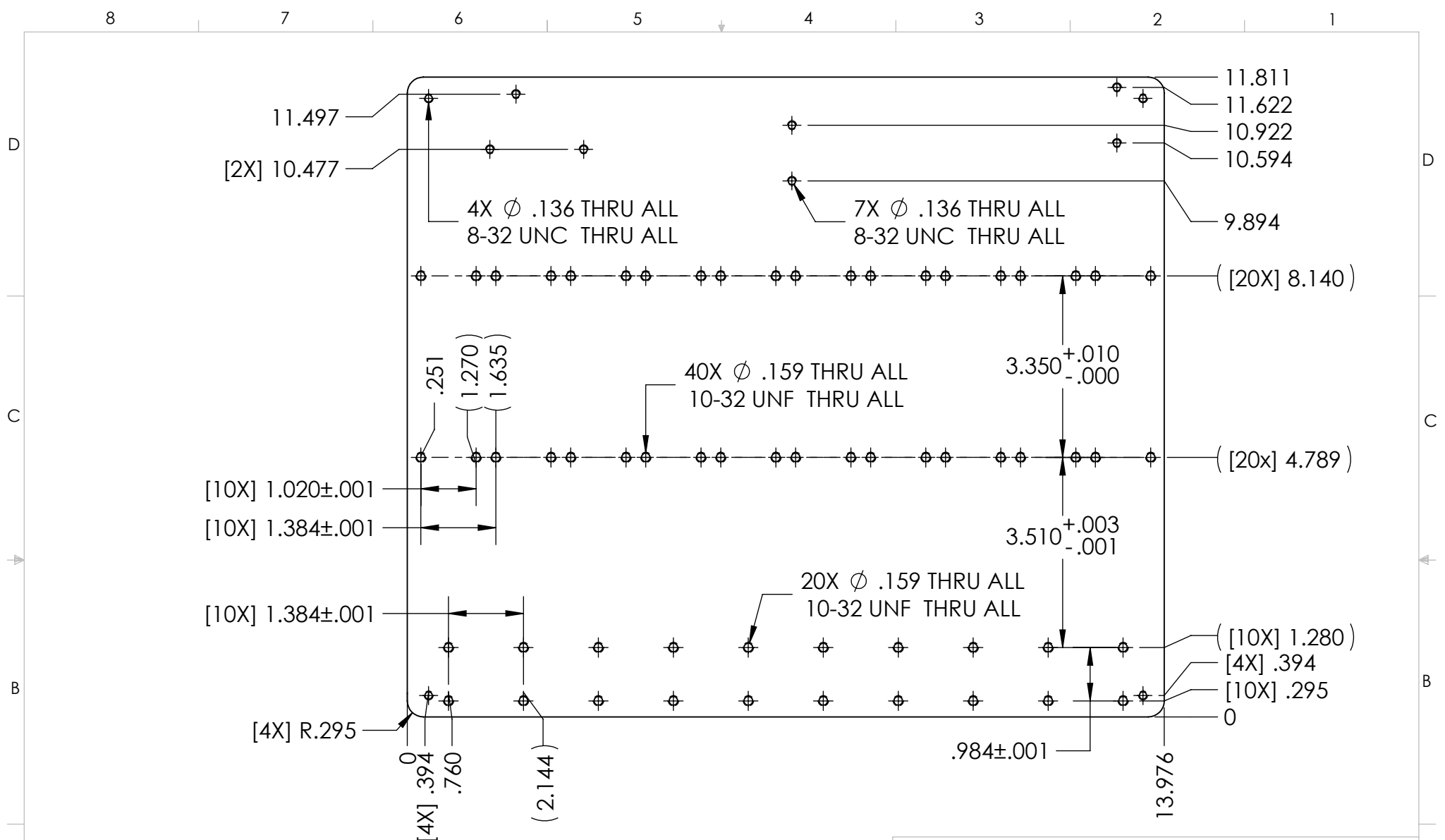

UNLESS OTHERWISE SPECIFIED:

1. INTERPRET DRAWING PER ASME Y14.5-2018. DIMENSIONS IN INCHES.
2. REMOVE ALL BURRS AND BREAK SHARP EDGES .005-.015.
3. ALL TAPPED HOLES TO BE 2B CLASS OF FIT, FULL THREAD THRU.
4. TAP DRILL DIAMETERS SHOWN ARE REFERENCE ONLY.
5. PART TO BE FREE OF CHIPS, CUTTING FLUID, AND CONTAMINATION PRIOR TO SHIPMENT.
6. MATERIAL: [ALUMINUM MIC-6] FINISH: [N/A]





CONFIDENTIAL			
TOLERANCES AND INSPECTION DIMENSIONS INCLUDED IN ASSOCIATED SPECIFICATION DOCUMENT		Q'APEL MEDICAL DO NOT REPRODUCE WITHOUT WRITTEN PERMISSION OF Q'APEL MEDICAL	
DIMS IN OVALS TO BE INSPECTED			
MATERIAL : ALUMINUM MIC-6		TITLE : FIXTURE BASEPLATE	
FINISH : N/A	DRAWN BY : RAVI NATARAJ	DOC NO:	REV:
CAD FILE :			2 OF 5

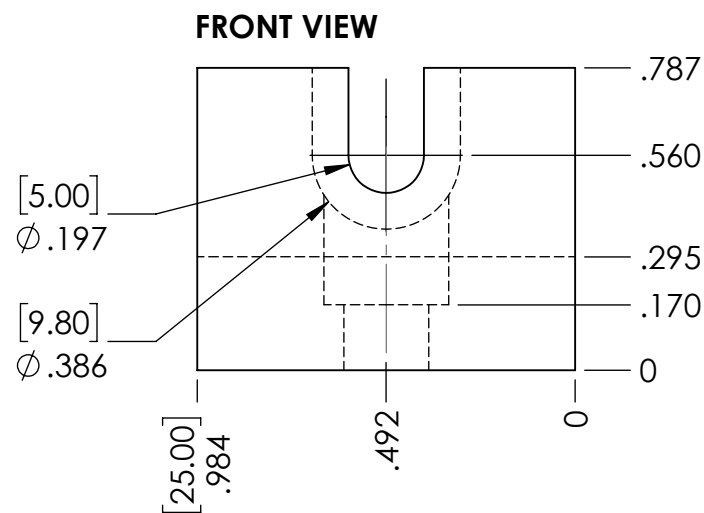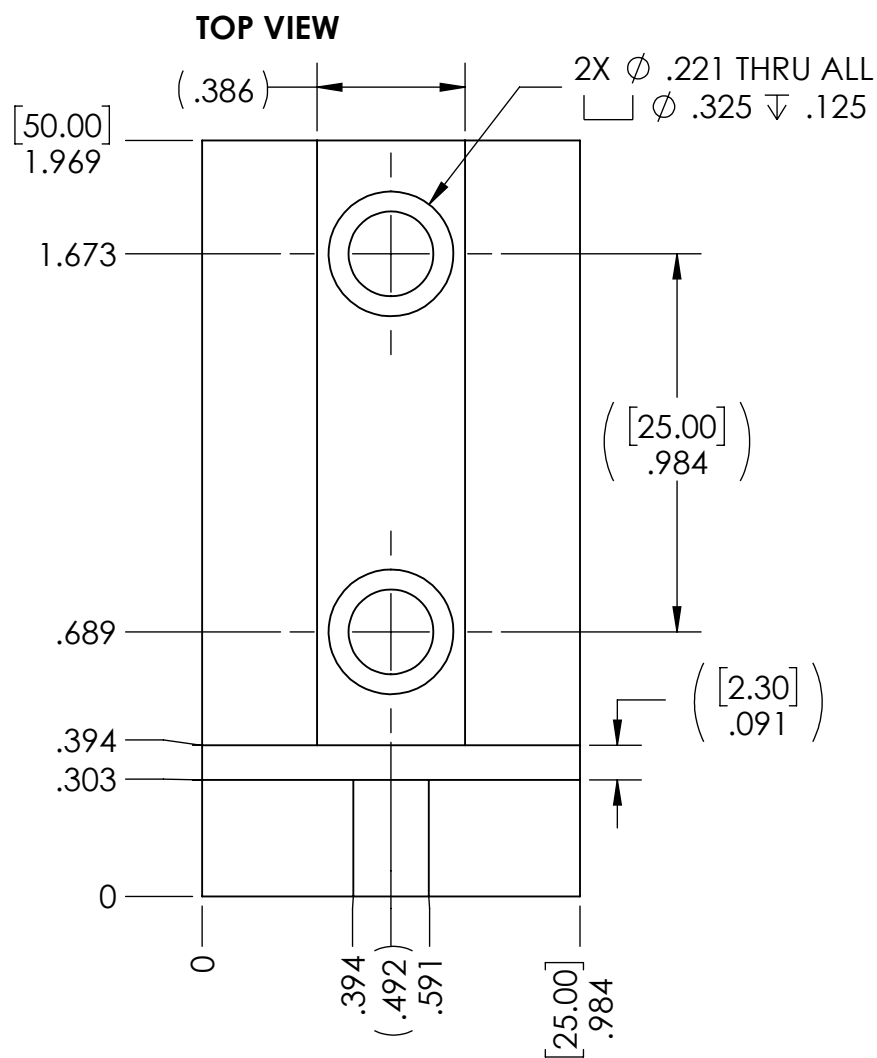

UNLESS OTHERWISE SPECIFIED:

1. DIMENSIONS IN INCHES. TOLERANCE $\pm .020$. NO INSPECTION REQUIRED.
2. 3D PRINTED, FDM. MATERIAL: [PLA/ABS/PETG].
3. REMOVE ALL SUPPORT MATERIAL.
4. DRILL $\varnothing .221$ HOLES AND $\varnothing .197$ BORE TO SIZE AFTER PRINTING IF REQUIRED

CONFIDENTIAL			
TOLERANCES AND INSPECTION DIMENSIONS INCLUDED IN ASSOCIATED SPECIFICATION DOCUMENT		Q'APEL MEDICAL DO NOT REPRODUCE WITHOUT WRITTEN PERMISSION OF Q'APEL MEDICAL	
DIMS IN OVALS TO BE INSPECTED			
MATERIAL : PLA - 3D-PRINTED		TITLE : 1ML-SYRINGE- HOLDER	
FINISH : N/A	DRAWN BY : RAVI NATARAJ	DOC NO:	REV:
CAD FILE :			3 OF 5

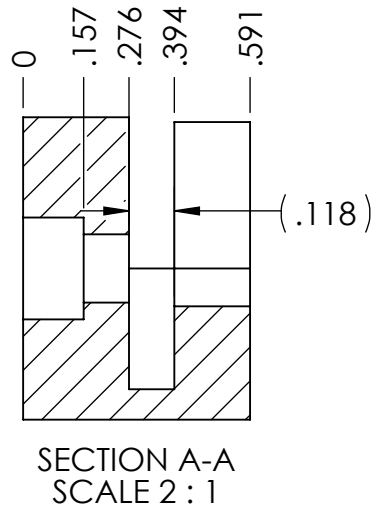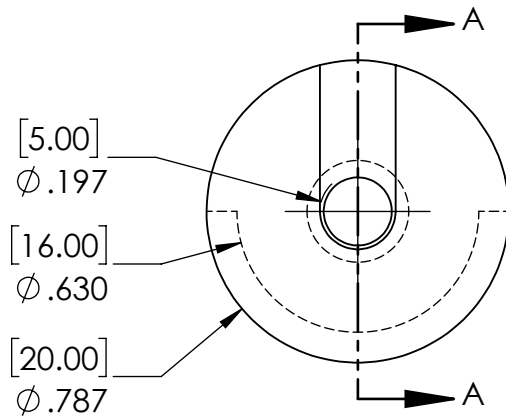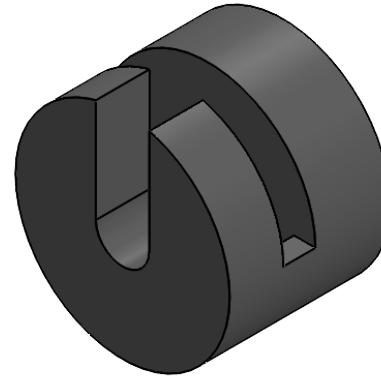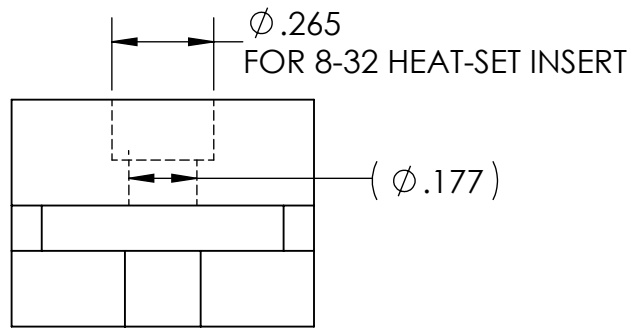


UNLESS OTHERWISE SPECIFIED:

1. DIMENSIONS IN INCHES. TOLERANCE $\pm .020$. NO INSPECTION REQUIRED.
2. 3D PRINTED, FDM. MATERIAL: [PLA/ABS/PETG].
3. PRINT WITH BORE AXIS VERTICAL. 4 PERIMETERS MIN.
4. REAM Ø.197 BORE AND DRILL Ø.265 INSERT HOLE TO SIZE AFTER PRINTING.
5. INSTALL 8-32 BRASS HEAT-SET INSERT PER MFR INSTRUCTIONS.
6. REMOVE ALL SUPPORT MATERIAL.

CONFIDENTIAL

TOLERANCES AND INSPECTION DIMENSIONS INCLUDED IN ASSOCIATED SPECIFICATION DOCUMENT		Q'APEL MEDICAL DO NOT REPRODUCE WITHOUT WRITTEN PERMISSION OF Q'APEL MEDICAL	
DIMS IN OVALS TO BE INSPECTED			
MATERIAL : PLA - 3D-PRINTED		TITLE : 1-ML-SYRINGE-PUSHER	
FINISH : N/A	DRAWN BY : RAVI NATARAJ	DOC NO:	REV:
CAD FILE :			4 OF 5

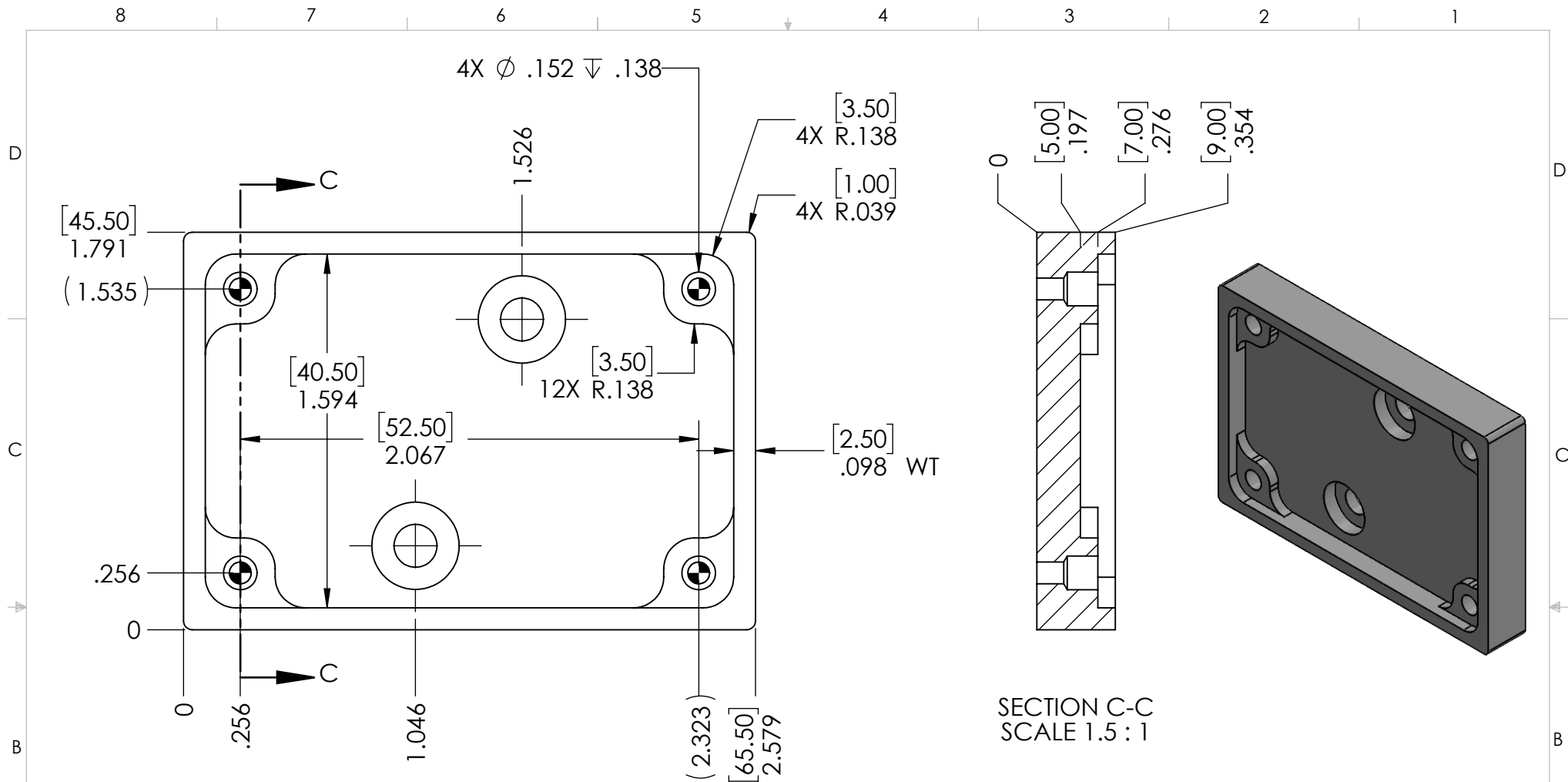

UNLESS OTHERWISE SPECIFIED:

1. DIMENSIONS IN INCHES [MILLIMETERS]. TOLERANCE $\pm$.020 [$\pm$ 0.5]. NO INSPECTION REQUIRED.
2. 3D PRINTED, FDM. MATERIAL: [PLA/ABS/PETG].
3. PRINT OPEN FACE DOWN ON BUILD PLATE. 4 PERIMETERS MIN.
4. INSTALL 4X M3 BRASS HEAT-SET INSERTS PER MFR INSTRUCTIONS.
5. PART MOUNTS DIRECTLY TO SOLENOID BODY. VERIFY FIT BEFORE FINAL ASSEMBLY.

CONFIDENTIAL			
TOLERANCES AND INSPECTION DIMENSIONS INCLUDED IN ASSOCIATED SPECIFICATION DOCUMENT		Q'APEL MEDICAL DO NOT REPRODUCE WITHOUT WRITTEN PERMISSION OF Q'APEL MEDICAL	
DIMS IN OVALS TO BE INSPECTED			
MATERIAL : PLA - 3D-PRINTED		TITLE : PCB-MOUNT ON 5/2 SOLENOID	
FINISH : N/A	DRAWN BY : RAVI NATARAJ	DOC NO:	REV:
CAD FILE :			5 OF 5